

ONYX[™]

Freedom Scientific

Swing-arm



ONYX Swing-arm at a Glance:

- Great agility - moves quickly between distance, document, and self views
- Turns your TV or VGA computer monitor into a video magnifier
- Find feature establishes image perspective before zooming and self views
- Magnification up to 100x (20" screen)

Three Distinct Views

This small, yet powerful, camera provides magnification up to 100x and is uniquely mounted on a swinging arm to provide you with three distinct views: distance, document, and self views. Easily connects to your TV or computer monitor for an affordable video magnification solution.

Putting You in Control

Change the magnification, brightness, contrast, and tracking features with a

unique and easy-to-use RF remote control. There is no need to touch the camera once it has been set, and you don't even need to point the remote control at the camera as with an infrared remote control.

The 11 3/8-inch workspace under the camera provides the room you need to work on objects with your hands and tools. Use the Find Feature to aim at a distant object before zooming in. Auto Focus eliminates the need for constant refocusing as the viewing field changes. Use Focus Lock to keep an object or document in focus when working or writing under the camera. The unique Position Sensor function retains



ONYX Swing-arm in self view

your settings in each of the three modes — distance view, document view, and self view — without need to readjust the camera as you move from one view to another and back again.

The many features of the ONYX Swing-arm include:

- Five selectable high-contrast modes - four modes when connected to a TV
- Brightness adjustments to eliminate glare and to enable viewing in low-light environments


- Masking and Reading Lines
- Find Feature
- Focus Lock
- Freeze Frame
- RF remote control
- Detachable camera
- Universal 110-240 volt power supply

Specifications

Weight:	5.7 lbs / 2.6 kg	Magnification at 113/8 in / 28.9 cm:
Camera-only weight:	14.8 oz / 6.7 kg	1.5x to 100x*
Working space under camera:	11 3/8 in / 28.9 cm	Magnification at 3 ft / 91.4 cm:
		25x*

*Based on 20-inch TV

Designed & Assembled in the USA 

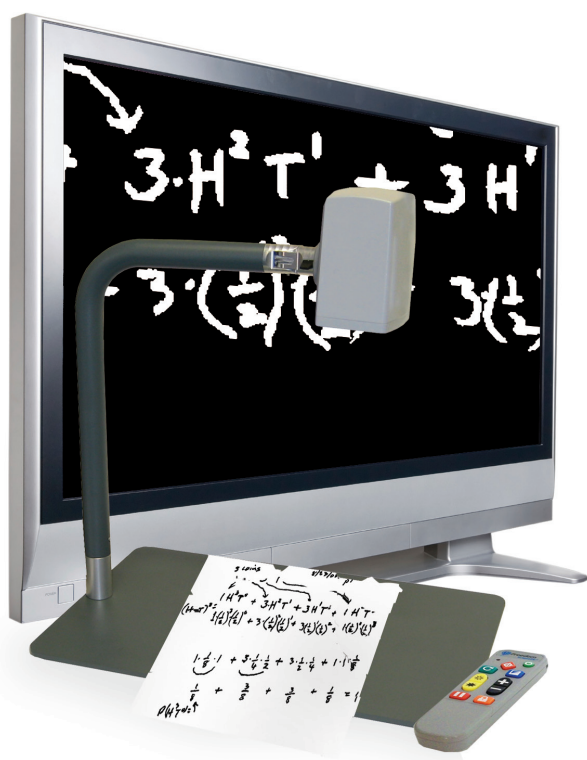

About Freedom Scientific

Freedom Scientific is the world's leading manufacturer of assistive technology products for those who are blind or vision impaired. For more information, e-mail info@freedomscientific.com